



TECHNISCHE
UNIVERSITÄT
WIEN
Vienna University of Technology

VIENNA UNIVERSITY OF TECHNOLOGY

FACULTY OF PHYSICS

LABORATORY III

Laboratory Report

Röntgen

Authors:

Raul Wagner

Martin Kronberger

Group 301

Supervisor:

Someone

conducted on:
21 May 2025

Measurement Setup and Preparations

Setup

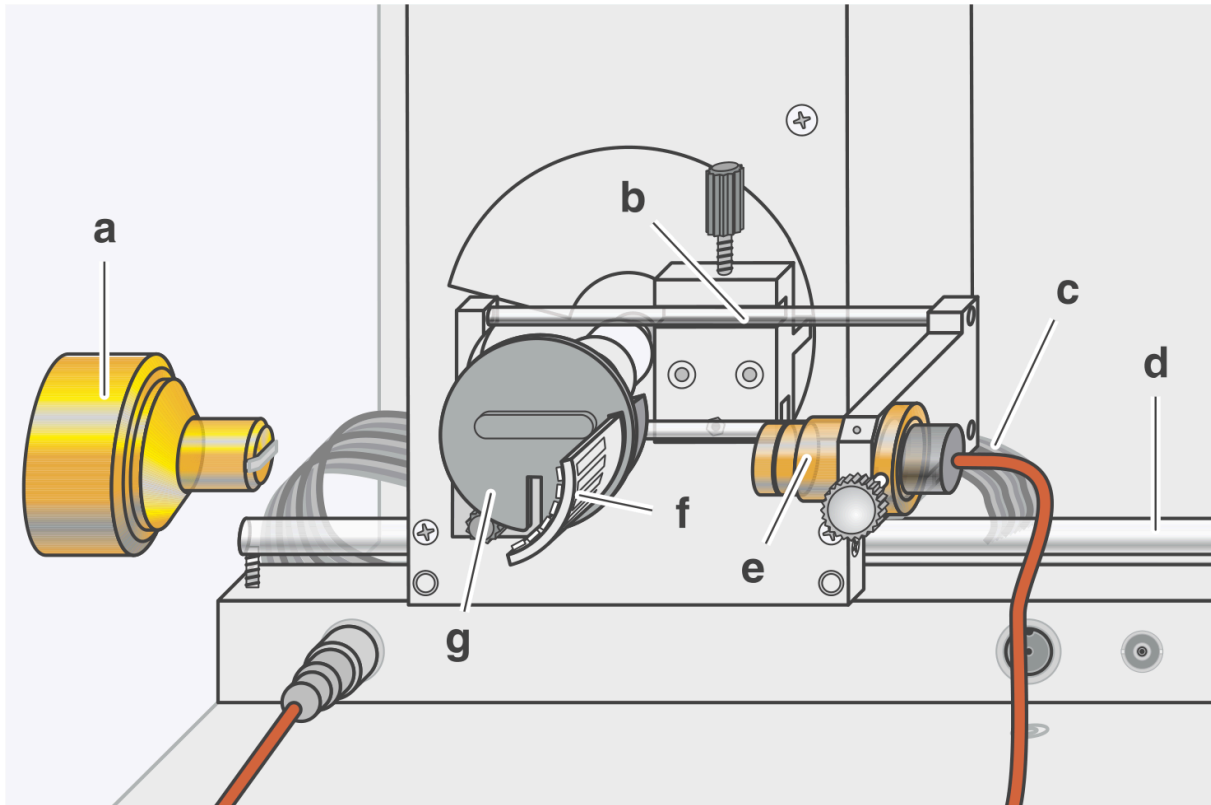


Figure 1: Measurement Setup with following components: (a) collimator mount, (b) sensor holder, (c) flat ribbon cable, (d) goniometer guide rods, (e) sensor mount, (f) insertion edge of absorber set I, and (g) goniometer target holder.

Preparations

- Carefully align the guide rod while inserting the collimator into the collimator mount (a).
- Secure the goniometer onto the guide rods (d) before connecting the flat ribbon cable (c) for control.
- After removing the protective cap, install the window counter tube into the sensor mount (e) and plug its cable into the GM-tube socket in the experimental area.
- Remove the goniometer's target holder (g) to lift off the target table.
- Slide the insertion edge of absorber set I (f) into the quarter-circle groove of the target holder until it clicks into place.
- Swap out the sensor holder with X-ray energy detector for the holder equipped with the window counter tube.
- Reinstall the target holder carrying absorber set I.
- Press the "Zero" button to set target and sensor to their null positions.
- Verify (and adjust if needed) the zero position of both the blank aperture in the absorber set and the sensor (see "Setting the measurement zero position" in the X-ray manual).
- Finally, slide the goniometer to position the collimator at 5 cm from the blank aperture, then slide the sensor holder (b) to set 5 cm between aperture and sensor slit. = Dependence of attenuation on absorber thickness

Experiment 1: Attenuation of X-Ray Radiation

Objectives

- Investigate the attenuation of X-ray intensity as a function of absorber thickness and material.
- Verify Lambert's law of exponential attenuation.
- Demonstrate the wavelength dependence of the attenuation coefficient.

Theory

When a narrow beam of X-rays of initial count rate R_0 passes through an absorber of thickness x , the transmitted count rate R satisfies

$$T = \frac{R}{R_0} \quad \text{where} \quad T(x) = e^{-\mu x} \quad \rightarrow \quad \ln T = -\mu x$$

where μ is the linear attenuation coefficient.

Setup

1. Mount the collimator and goniometer on the X-ray tube as shown in Fig. 1.
2. Insert the Geiger-Müller detector in the sensor arm and connect via "GM Tube".
3. Align the target (absorber holder) and detector so that the slit-to-target and target-to-detector distances are each ≈ 5 cm.
4. Zero-position both arms with the "Zero" button.

Attenuation vs. Absorber Thickness

Without zirconium filter

1. Set tube voltage $U = 21$ kV, emission current $I = 0.05$ mA, measurement time $\Delta t = 100$ s.
2. Set absorber angles corresponding to the corresponding thicknesses, press **Scan**, wait Δt , then read count rate R via **Replay**.
3. Record in Table 1.

| d / mm | R / s ⁻¹ |
|--------|---------------------|
| 0 | 1618 |
| 0.5 | 787.4 |
| 1 | 403.5 |
| 1.5 | 226.4 |
| 2 | 49.1 |
| 2.5 | 30.55 |
| 3 | 16.11 |

Table 1: Some Caption

With zirconium filter

1. Mount Zr filter, set $I = 0.15$ mA, $\Delta t = 200$ s.
2. Repeat step 1.1 at the same angles.
3. Record in Table 2.

| d / mm | R / s ⁻¹ |
|--------|---------------------|
| 0 | 775.1 |
| 0.5 | 337 |
| 1 | 149.8 |
| 1.5 | 72.1 |
| 2 | 33.85 |
| 2.5 | 18.6 |
| 3 | 7.85 |

Table 2: Some caption

Data Analysis

- Compute transmission: $T(d) = \frac{R(d)}{R(0)}$.
- Plot $T(d)$ vs d and $\ln T(d)$ vs d .
- Fit $\ln T(d) = -\mu d$ to extract μ for both unfiltered and filtered cases.

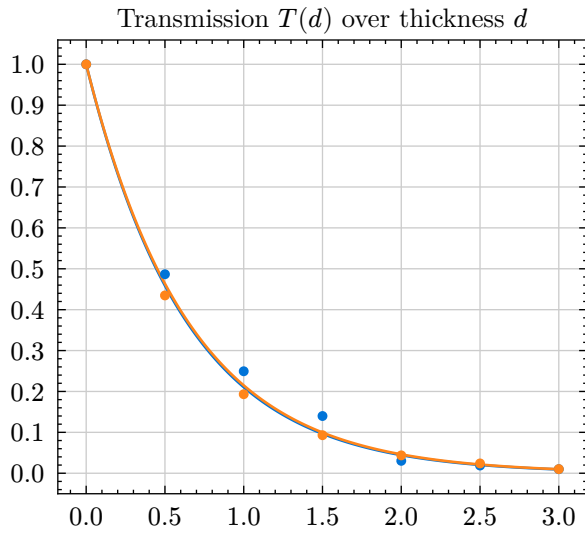


Figure 2:

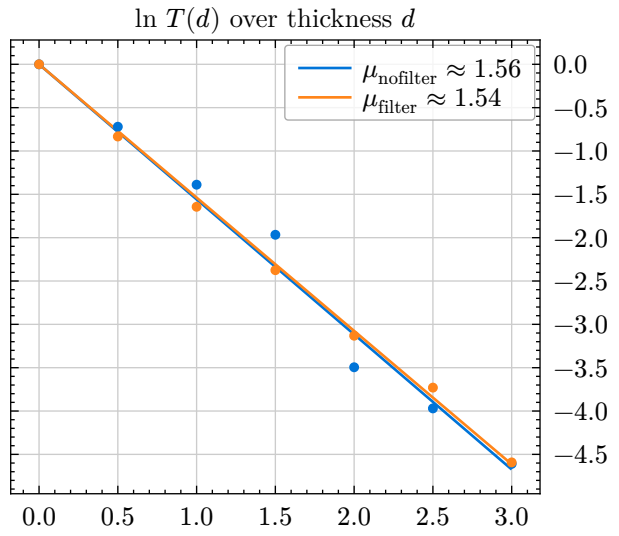


Figure 3:

Dependence of attenuation on the absorber material**Without zirconium filter**

| Absorber | Z | I / mA | Δt / s | R / s ⁻¹ | T | μ / cm ⁻¹ |
|----------|----|--------|----------------|---------------------|-------|--------------------------|
| leer | 0 | 0.02 | 30 | 89550 | 1 | 0 |
| C | 6 | 0.02 | 30 | 85050 | 0.95 | 17 |
| Al | 13 | 0.02 | 30 | 54150 | 0.605 | 168 |
| Fe | 26 | 1 | 300 | 7355 | 0.082 | 833 |
| Cu | 29 | 1 | 300 | 15.55 | 0 | 2886 |
| Zr | 40 | 1 | 300 | 181.8 | 0.002 | 2067 |
| Ag | 47 | 1 | 300 | 56.65 | 0.001 | 2455 |

Table 3: Some caption

With a zirconium filter

| Absorber | Z | I / mA | Δt / s | R / s ⁻¹ | T | μ / cm ⁻¹ |
|----------|----|--------|----------------|---------------------|-------|--------------------------|
| leer | 0 | 0.02 | 30 | 36915 | 1 | 0 |
| C | 6 | 0.02 | 30 | 34730 | 0.941 | 20 |
| Al | 13 | 0.02 | 30 | 19545 | 0.529 | 212 |
| Fe | 26 | 1 | 300 | 2585 | 0.07 | 886 |
| Cu | 29 | 1 | 300 | 5.7 | 0 | 2925 |
| Zr | 40 | 1 | 300 | 107.3 | 0.003 | 1947 |
| Ag | 47 | 1 | 300 | 12.05 | 0 | 2676 |

Table 4: Some caption

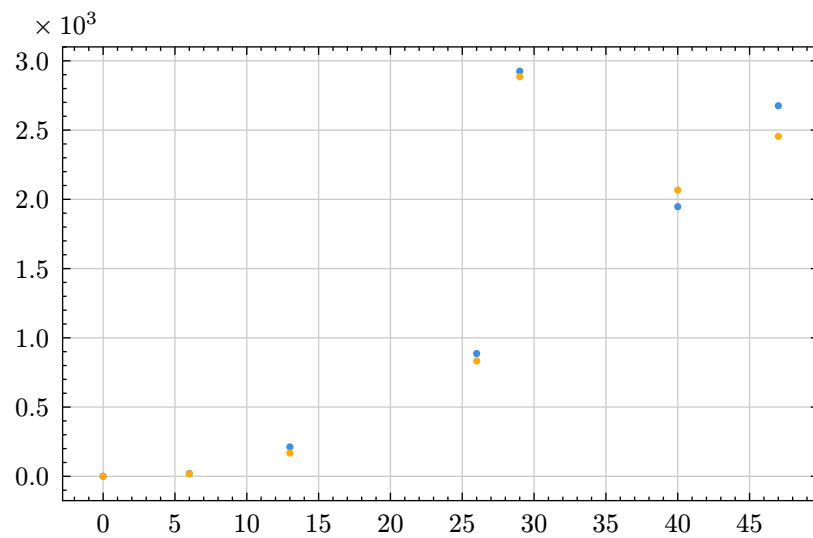


Figure 4: Some caption

Measurement of the Zeroeffect

Experiment 2: Bragg Reflection

Objectives

- Investigate the Bragg reflection of Mo K-characteristic X-rays on a NaCl single crystal.
- Determine the wavelengths of the $K\alpha$ and $K\beta$ lines up to third order diffraction.
- Confirm Bragg's law and the wave nature of X-radiation.

Theory

When X-rays hit parallel crystal planes spaced by distance d , constructive interference occurs at angles θ satisfying Bragg's law:

$$n \sin \theta = n \frac{\lambda}{2d}$$

where n is the diffraction order and λ is the wavelength.

Apparatus

- X-ray tube with collimator mount and goniometer guide rods.
- Geiger-Müller detector in the sensor mount.
- NaCl single crystal fixed on the crystal stage.
- Distances: collimator-crystal ≈ 5 cm; crystal-detector ≈ 6 cm.

Procedure

1. Connect PC via USB, start the "X-ray Device" software, select automatic scan mode.
2. Set parameters: tube voltage $U = 35$ kV, current $I = 1.00$ mA, measurement time $\Delta t = 10$ s, angle step $\Delta\beta = 0.1^\circ$; press COUPLED; set scan range $\beta_{\min} = 2.5^\circ$, $\beta_{\max} = 30^\circ$.
3. Start scan to record the spectrum; save data with F2.
4. Identify peak angles θ for $K\alpha$ and $K\beta$ lines at orders $n = 1, 2, 3$ and record in tables.
5. Calculate wavelengths via $\lambda = 2d \sin \theta$ using $d = 282.01$ pm.

Results

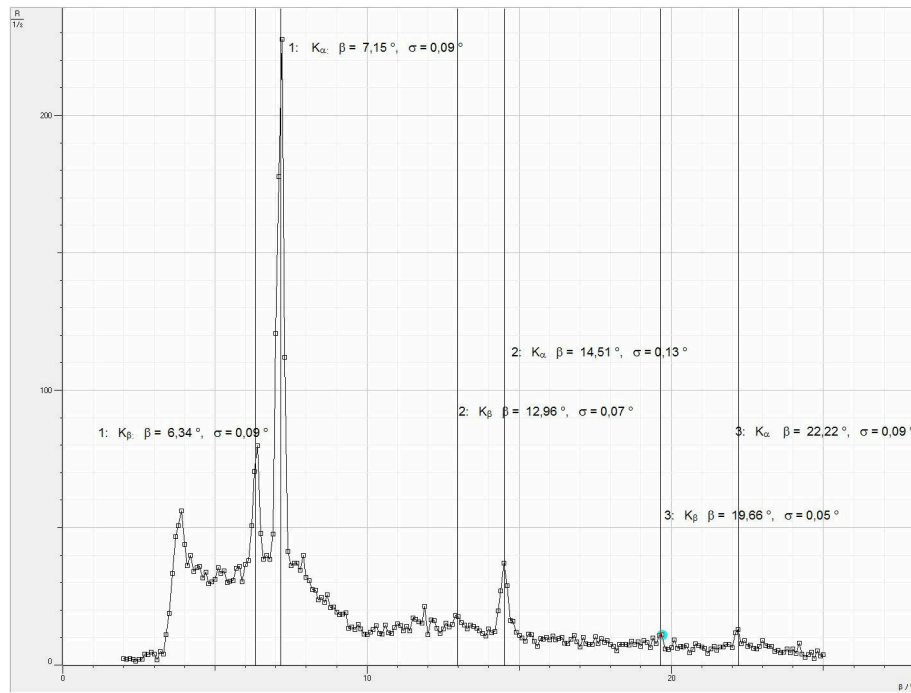


Figure 5: Some Caption

Conclusions

Measured wavelengths agree with literature and validate Bragg's law and the wave character of X-rays.

Experiment 3: Duane–Hunt Law and Planck’s Constant

Objectives

- Determine the cutoff wavelength λ_{\min} of the Bremsstrahlung continuum as a function of tube voltage U .
- Verify the Duane–Hunt relation $\lambda_{\min} = \frac{hc}{eU}$.
- Extract Planck’s constant h from the slope of λ_{\min} vs $\frac{1}{U}$.

Theory

Complete conversion of electron kinetic energy into photon energy gives:

$$\lambda_{\min} = \frac{hc}{eU}$$

where e is the elementary charge and c the speed of light.

Apparatus

- Same goniometer and NaCl crystal setup as in Experiment 2.
- “X-ray Device” software with Planck-mode register.
- Geiger–Müller detector.

Procedure

1. For tube voltages $U = 22, 24, 26, 28, 30, 32, 34, 35$ kV at $I = 1.00$ mA, set measurement time and angle range as in Table 1.
2. Perform automatic scans; save each spectrum.
3. In Planck mode, determine λ_{\min} for each U .
4. Plot λ_{\min} vs. $\frac{1}{U}$ and fit a line through the origin; extract slope A .

Results

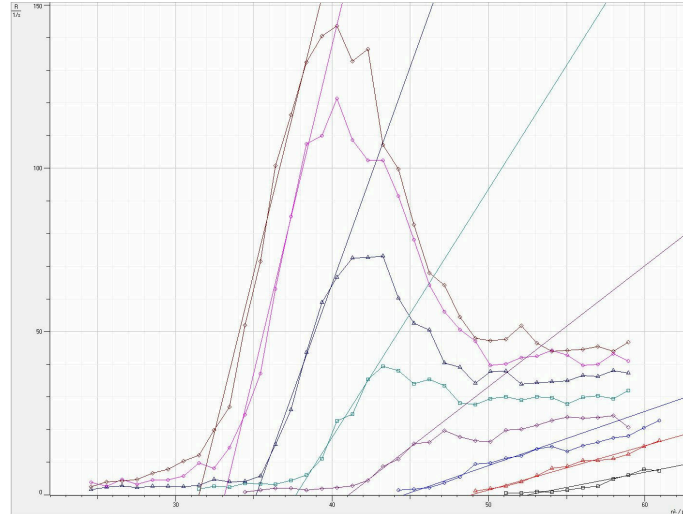


Figure 6: Some Caption

Laboratory Work III - Röntgen

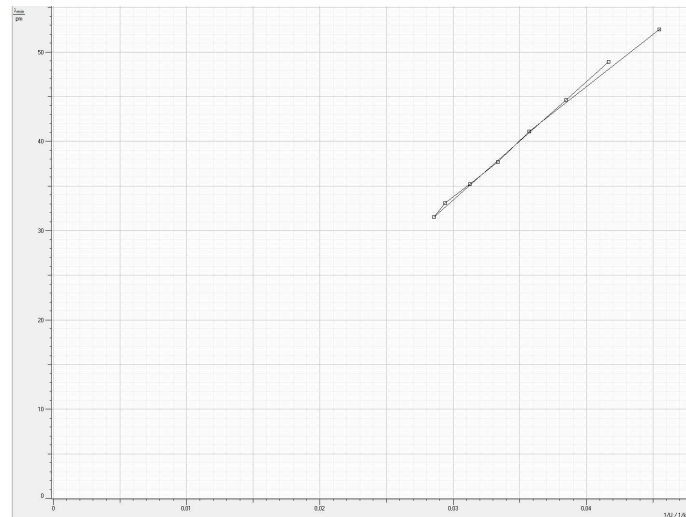


Figure 7: Some Caption

$$A = 1137 \text{ pm kV}$$

Conclusions

The Duane–Hunt law is confirmed, and the measured Planck constant agrees closely with the literature value.

Planks konstant as a wavelenth-voltage factor is about 1240 pm kV